


NASG Application: Facility Level


Always wear gloves when putting the NASG on a bleeding woman.




Place the NASG under the woman

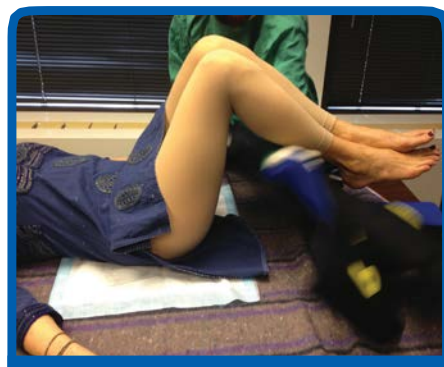

If she can help, first have her raise her feet

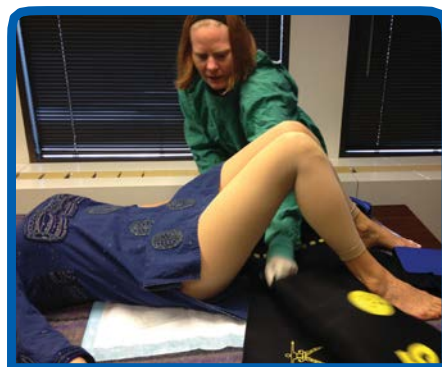


Then place her feet on the bed (not the NASG)

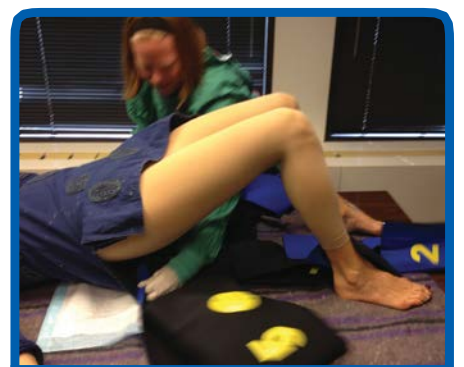

Have her raise her hips while you slide the NASG under her



Place the ball in segment #5 over her navel for fit, but do not close segments #5 and #6 until the end of the application

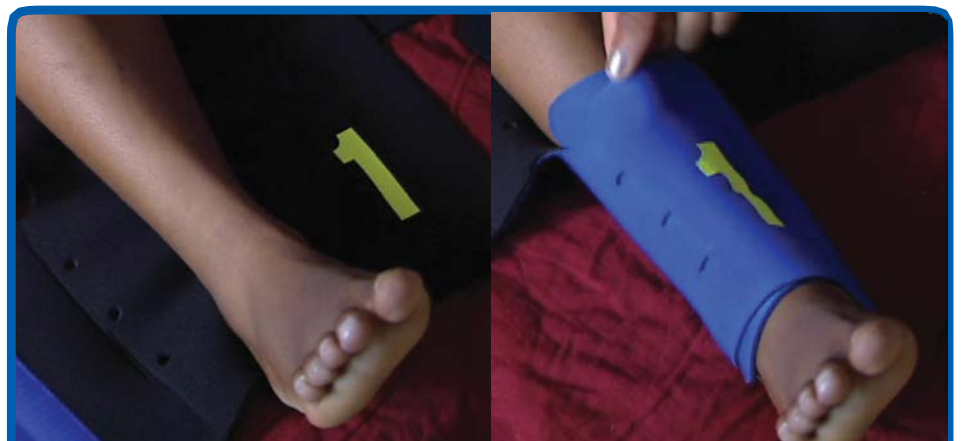

Start at the ankles, close segment pair #1 (or #2 for short women) around her ankles. Check for tightness by pulling up on each segment with your finger, let go, it should make a sharp




Close segment pair #2 as tightly as possible; snap each segment



Close segment pair #3 as tightly as possible; snap each segment



Close segment #4 around the woman's pubic bone and pelvis

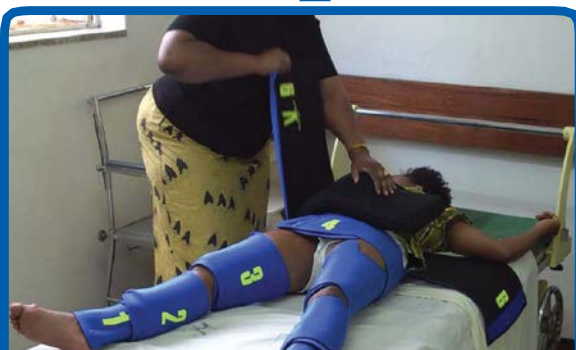

Put segment #5, with the ball, over the woman's navel



Close segments #5 and #6 as tightly as possible



Ask the woman if she can breathe normally

Continue emergency care of woman with hypovolemic shock and obstetric haemorrhage